

# A novel airborne nanoparticle size distribution spectrometer

H. Gnewuch, B. Gorbunov, R. Muir, R. Gilham\* and P. Quincey\*

Naneum Limited, CEH, University of Kent, Canterbury, Kent CT2 7NJ, UK

\* National Physical Laboratory, Hampton Road, Teddington, Middlesex, TW11 0LW, UK

A novel portable instrument SAC to record nanoparticle size distributions has recently been developed. It is based upon the principle of diffusion selection of particles drawn through a semi-penetrable net at various flow rates. The SAC instrument contains a size selective unit, condensation unit and an optical detector of particles. The instrument does not require charging particles, therefore, there are no uncertainties associated with the charge transfer process. The instrument enables nanoparticle size distributions to be recorded *in-situ* on-line in the range from 5 nm to 500 nm ([www.naneum.com](http://www.naneum.com)).

## Description

The principle of operation is as follows: Aerosol particles enter the size selective unit through the inlet. In the unit some particles are captured by the diffusion element. The number concentration of aerosol particles exiting the outlet of the size selective unit is measured subsequently by a particle counter and recorded by a PC. The number of particles captured by the size selective unit is influenced by the flow rate through it and the size distribution. The penetration of aerosol particles through net or screen type diffusion element is described in many publications, e.g. Yeh Hsu-Chi and Cheng Yung-Sung (1980). The variable flow rate is maintained by a pump, controlled by the PC and recorded by a mass flow sensor. Therefore, the instrument provides a number of measurements of particle concentrations at various flow rates.

A novel highly effective technique for aerosol particle size distribution reconstruction has been developed. It enables complicated size distributions to be obtained from measurements with the instrument developed.

## Preliminary tests results

The SAC instrument has been tested using atmospheric aerosols including products of combustion and some engineered nanoparticles. The engineered nanoparticles were produced by different generators: (i) electrosprayed sucrose, (ii) tungsten obtained with a hot-wire generator, and (iii) NaCl using an atomizer aerosol generator. Aerosols obtained were characterised by the novel instrument and reference methods. The SAC instrument was compared in a series of comparative measurements with TSI hand-held CPC and TSI SMPS.

For the experiments undertaken, the size distributions varied having a median diameter between ~15 and ~80 nm and a geometric standard deviation between ~1.4 and ~2.3. The comparison between the TSI SMPS measurements and the corresponding measurements obtained with the novel instrument showed a ratio between the median diameters to be  $\sim 1.2 \pm 0.2$  and the ratio between the geometric standard deviations to be varying between ~0.85 and ~1.2, see Figure 1. Meanwhile, the total number concentration measured with our instrument was at least ~80%

compared to the measurement with the TSI hand-held CPC. Therefore, the developed portable instrument can be used to characterise nanoparticle size distributions.

## References

Yeh Hsu-Chi and Cheng Yung-Sung (1980) Theory of a Screen-Type Diffusion Battery, *J. Aerosol Sci.*, **11**, pp. 313-320.

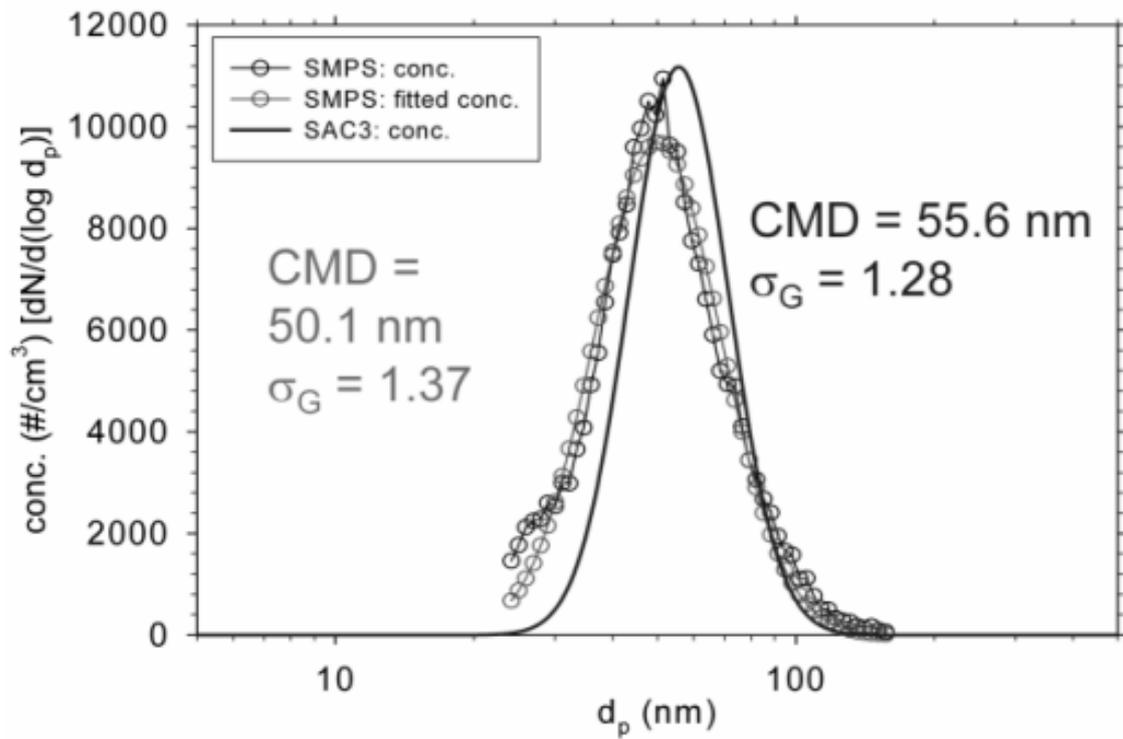


Figure 1. A typical aerosol size distribution (Tungsten) obtained with SAC NANEUM and SMPS TSI.

# A novel airborne nanoparticle size distribution spectrometer

H. Gnewuch, B. Gorbunov, R. Muir, R. Gilham\* and P. Quincey\*

Naneum Limited, CEH, University of Kent, Canterbury, Kent CT2 7NJ, UK

\* National Physical Laboratory, Hampton Road, Teddington, Middlesex,  
TW11 0LW, UK

# The technology of the size spectrometer

A novel portable instrument to record aerosol size distributions SAC is based upon the principle of **diffusion selection** of particles drawn through a semi-penetrable net at **various flow rates**

# Specification

- The size range - 2 nm to 500 nm equivalent diameter
- Measurement time – 20 s to 3 min
- Flow rates – 0.05 to 3 l/min
- Dimensions ( LWH , cm) - 30 x 20 x 23
- Weight - 5 k g
- Power – 5 NiMH batteries





# Main elements

The SAC instrument contains:

- a size selective unit,
- condensation unit and
- an optical detector of particles

# Schematic of SAC aerosol size spectrometer

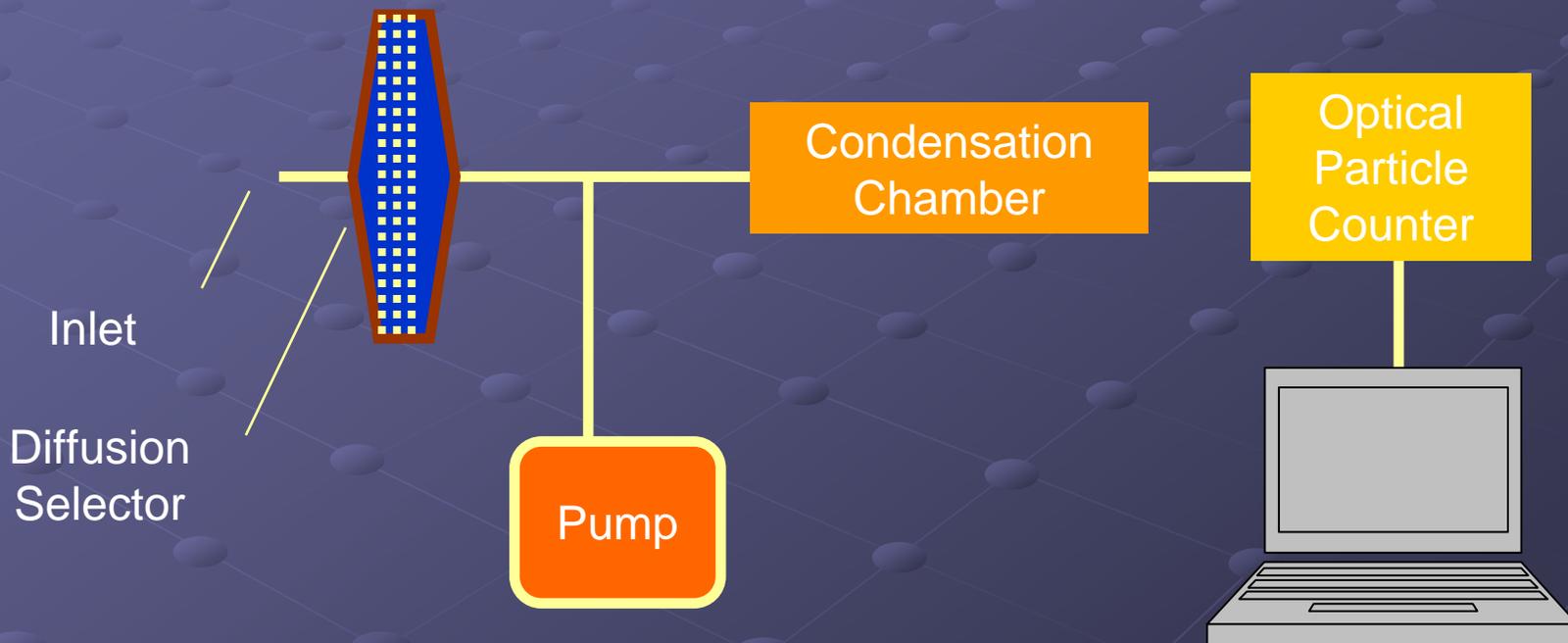


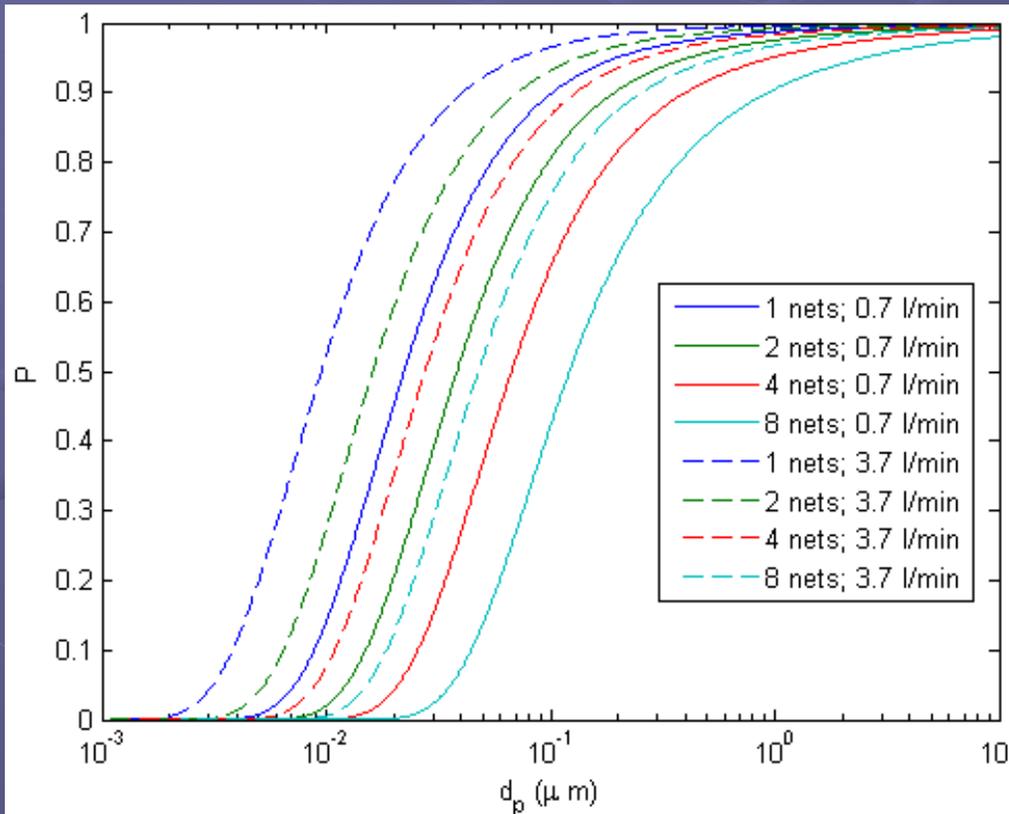


## Theory of particle deposition

Penetration  $P$  versus particle diameter  $d_p$  for various numbers of screens (nets) and flow rates (only diffusion as underlying deposition principle considered)

# The size selection is based upon diffusion

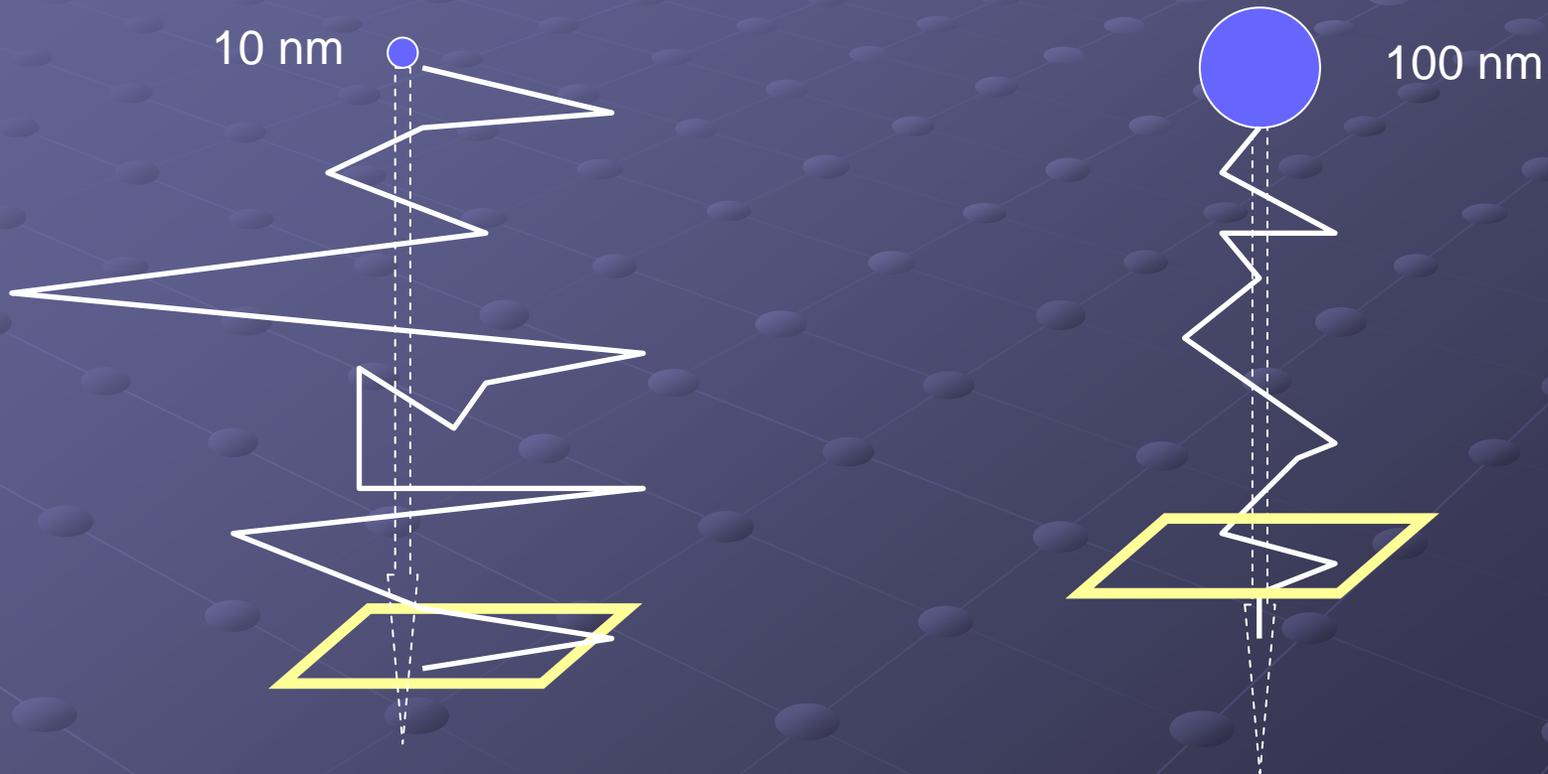


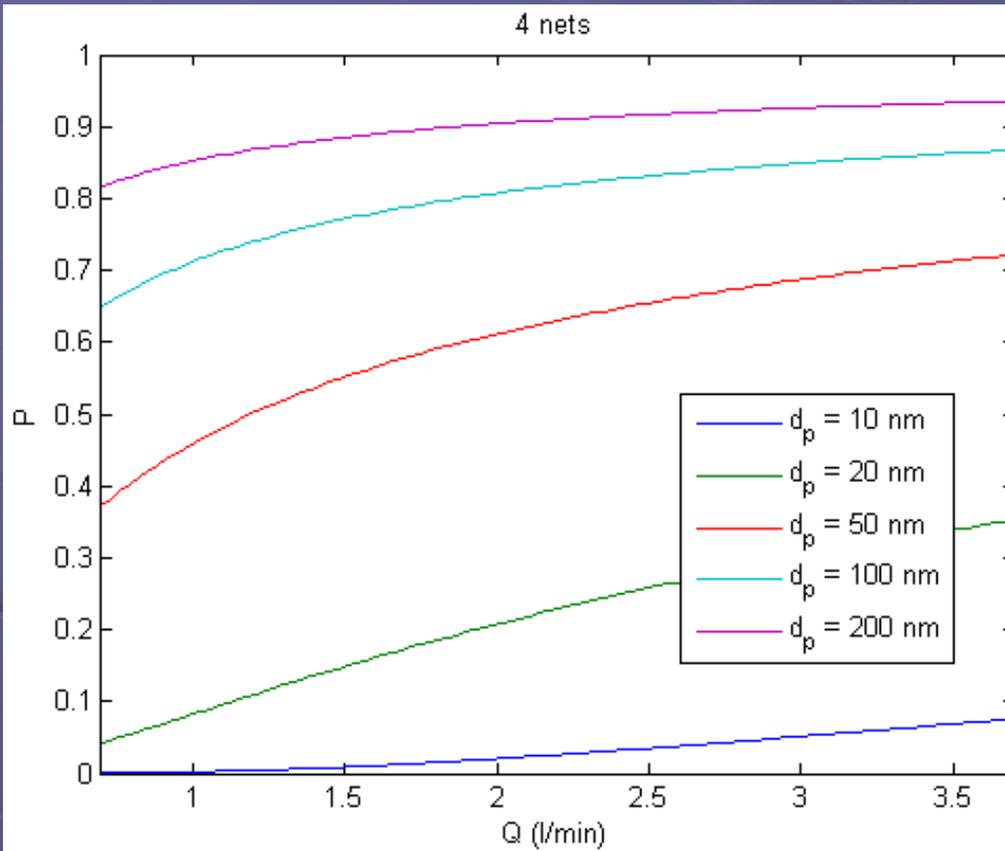


Theory of particle deposition for SAC aerosol particle size spectrometer

Penetration P versus flow rate Q for different particle diameters


# Schematic of SAC aerosol size spectrometer

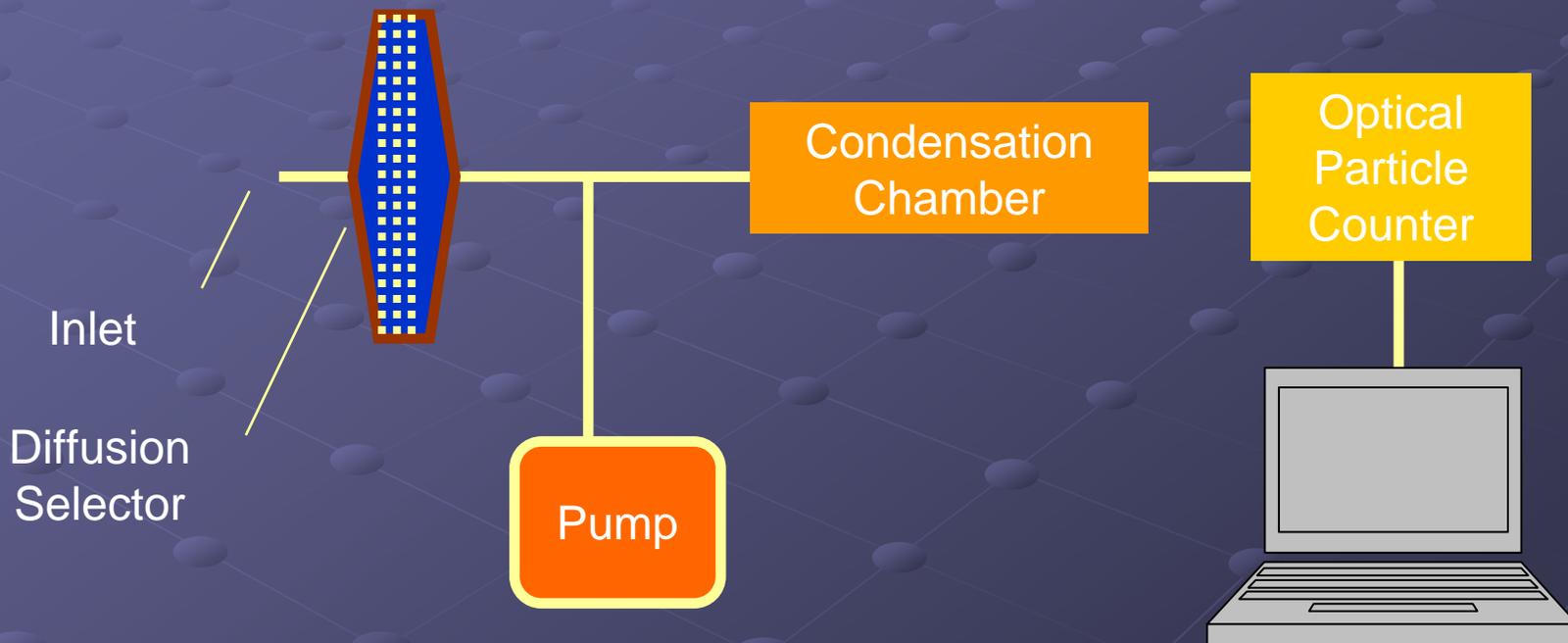


# Comparison: SAC vs. SMPS

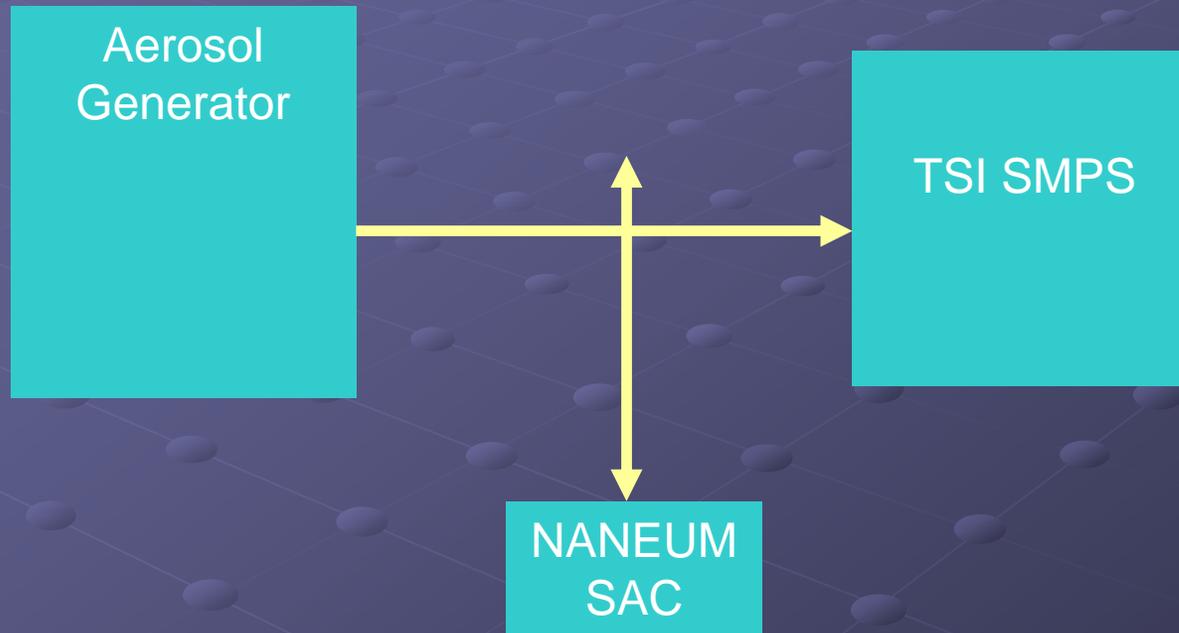


# Measurement protocol at NPL :

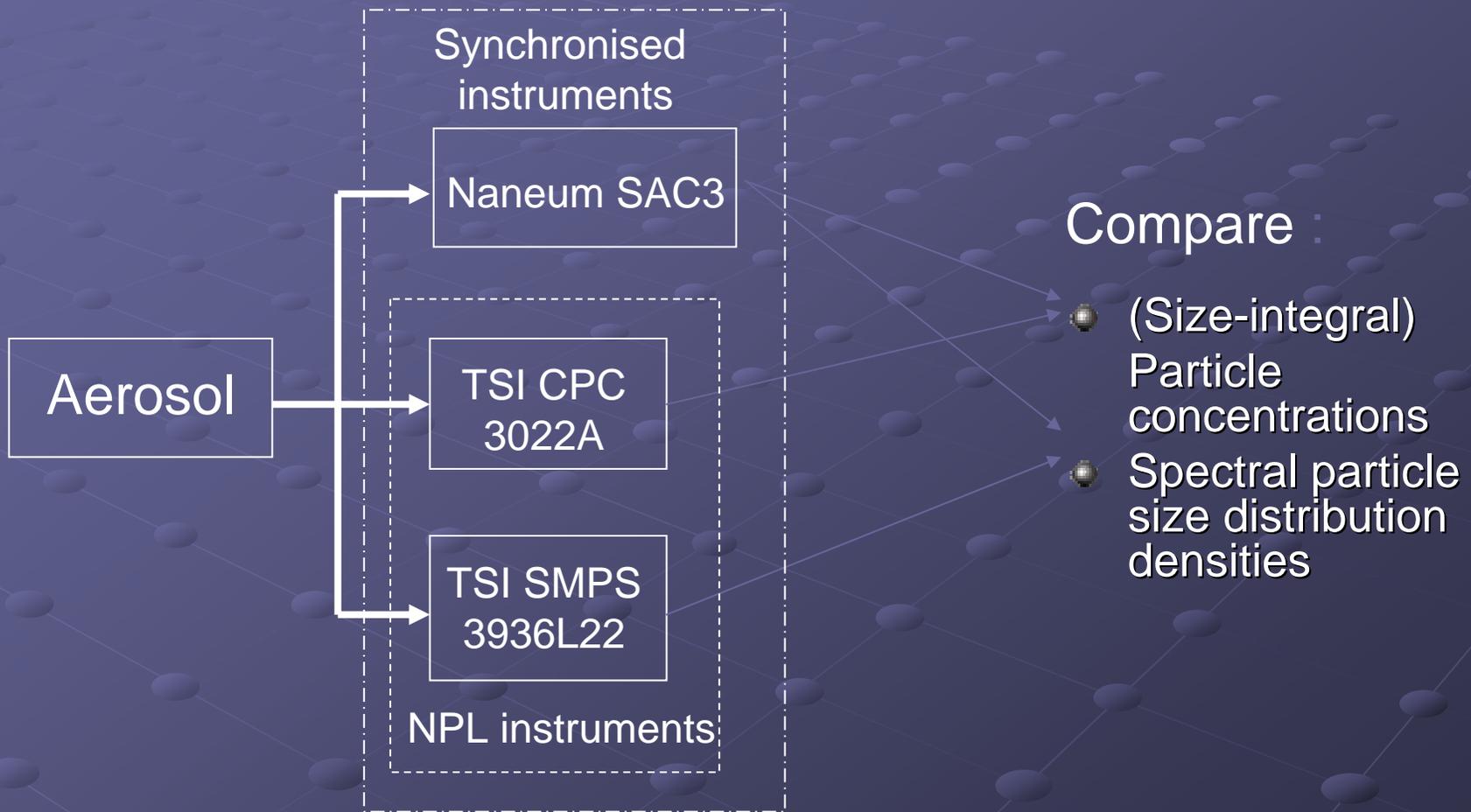


# Comparison with SMPS TSI (NPL) Sucrose

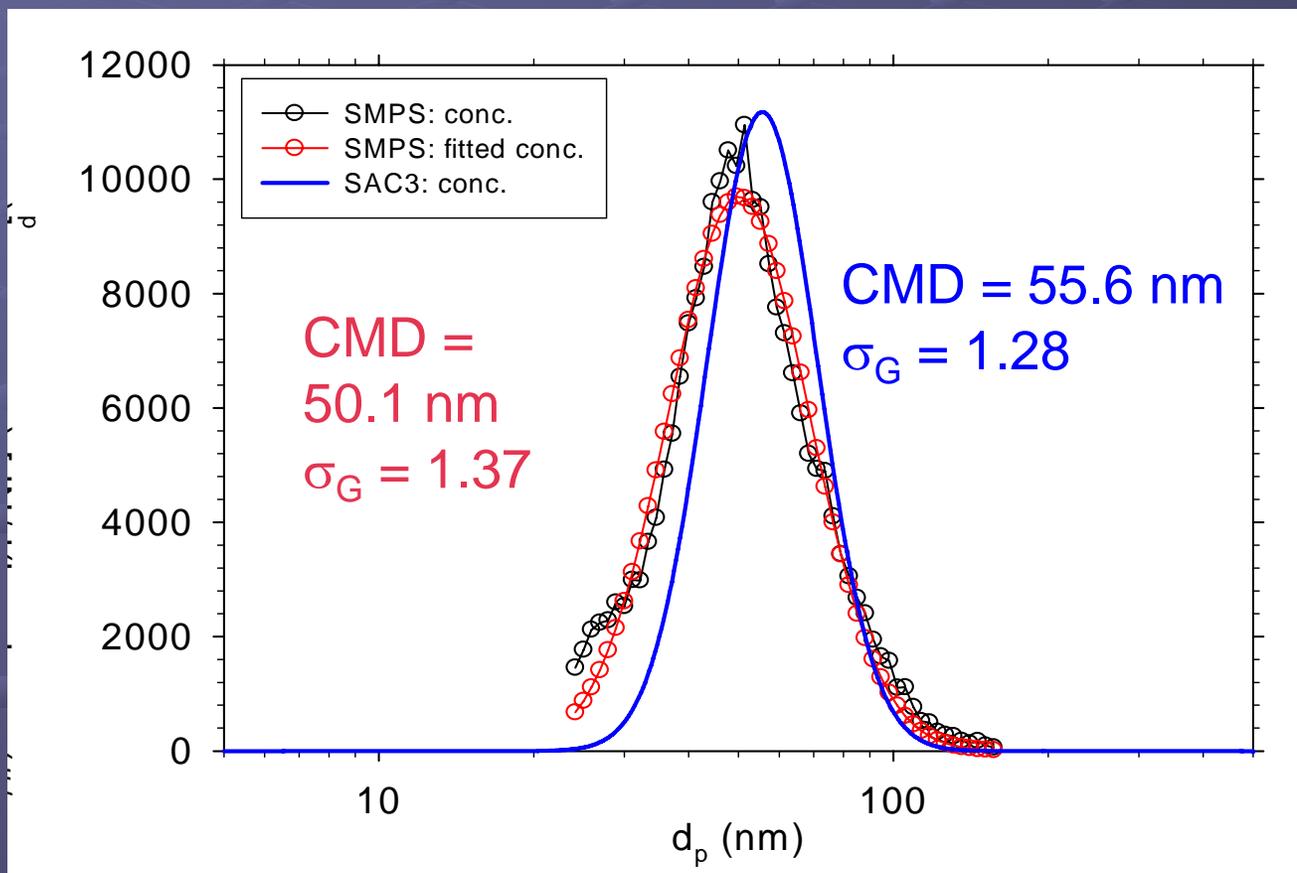

# Sodium Chloride aerosol generated by an atomiser

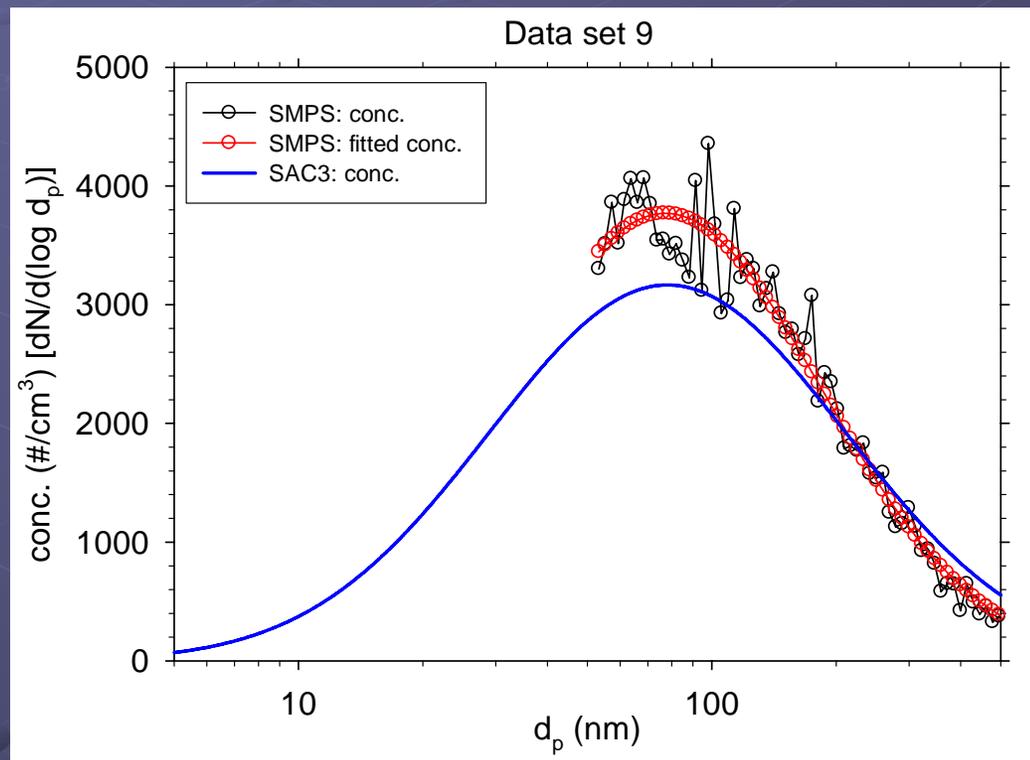

# Tungsten aerosol

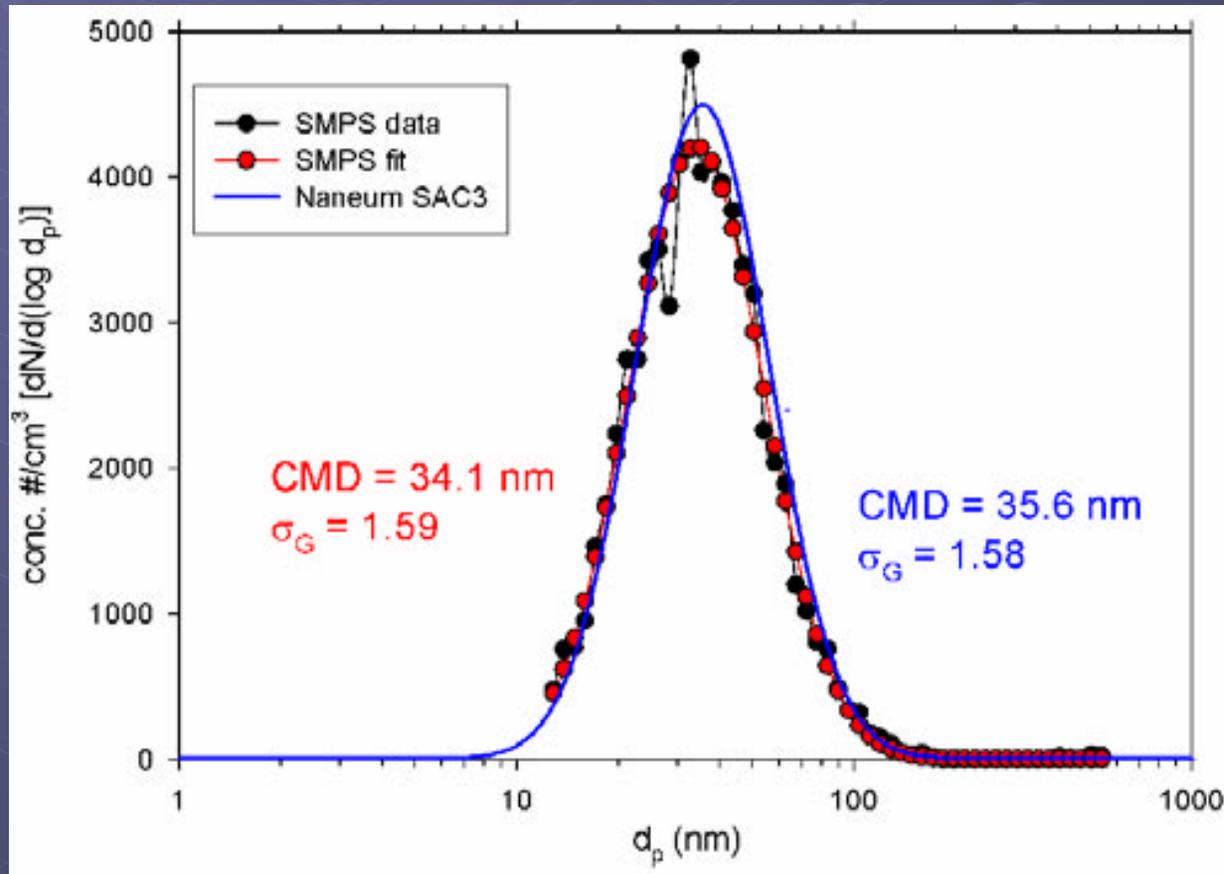

## Preliminary tests results

- The SAC instrument has been tested using atmospheric aerosols including products of combustion and some engineered nanoparticles (median diameter between ~15 and ~80 nm and a geometric standard deviation between ~1.4 and ~2.3).
- The engineered nanoparticles were produced by different generators: (i) electrosprayed sucrose, (ii) tungsten obtained with a hot-wire generator, and (iii) NaCl using an atomizer aerosol generator.
- The SAC instrument and TSI SMPS were in a good agreement of quantifying aerosol number particle size distributions

Naneum SAC 3.1

File Run Help

Set Name:

Number of Screens:


Number of Steps:

Step Duration (sec.):

Comment

Solutions

Atmospheric aerosol

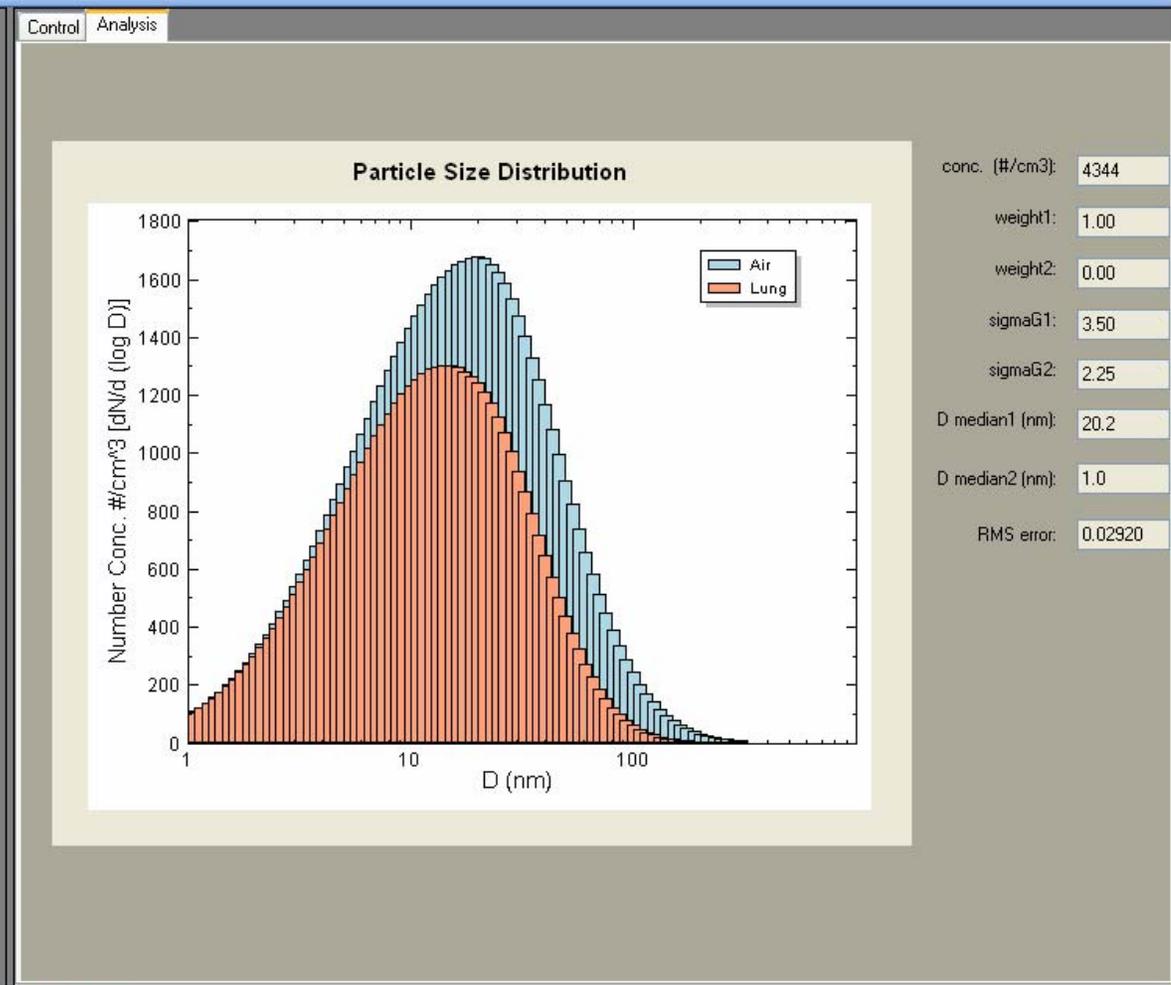

ConnectivityStatus ReadyStatus SteppingStatus AnalysisStatus

**Naneum SAC 3.1**

File Run Help

Set Name:

Number of Screens:

Number of Steps:

Step Duration (sec.):

Comment

Solutions

Control Analysis

### Particle Size Distribution

Number Conc. #/cm<sup>3</sup> [dN/d (log D)]

D (nm)

(1.273, 856.8)

Legend: Air, Lung

conc. (#/cm<sup>3</sup>):

weight1:

weight2:

sigmaG1:

sigmaG2:

D median1 (nm):


D median2 (nm):

RMS error:

ConnectivityStatus ReadyStatus SteepingStatus AnalysisStatus

start KINGSTON (E:) sac1.bmp - Paint EN 11:12

### Particle Size Distribution

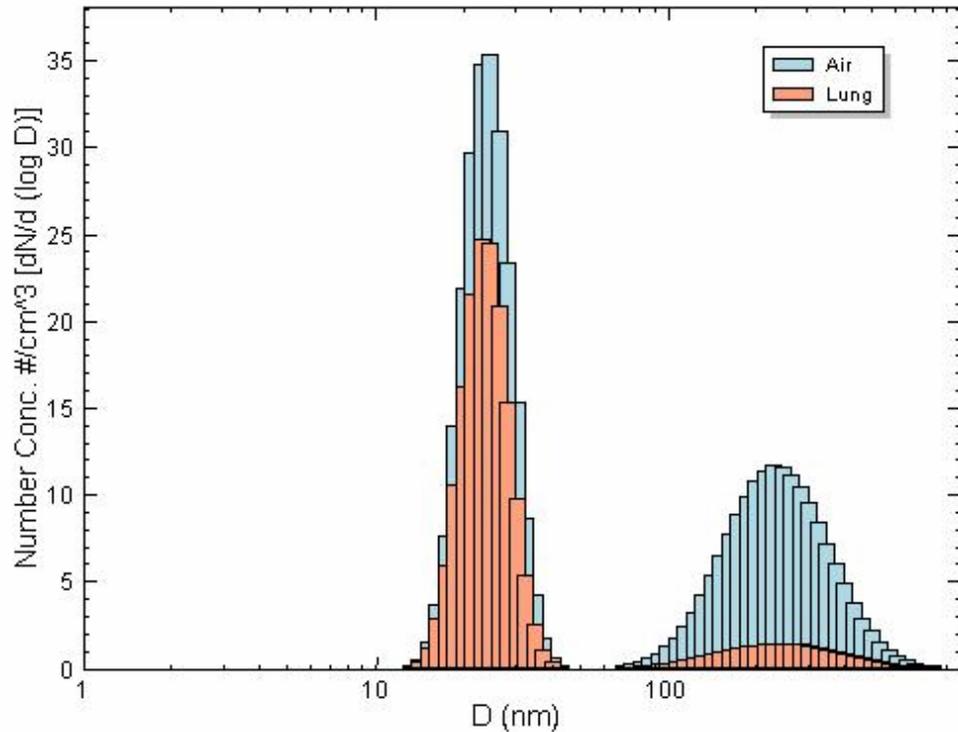

An example of bi-modal size distribution measured with SAC Aerosol Size Spectrometer

## Conclusions

1. The comparison showed a ratio between the median diameters to be  $\sim 1.2 \pm 0.2$  and the ratio between the geometric standard deviations to be varying between  $\sim 0.85$  and  $\sim 1.2$
2. The instrument developed can be used to quantify complex, bimodal aerosol particle size distributions *in situ* in real time
3. The instrument is robust, reliable and easy to use
4. There is no ionisation stage in measuring size distributions and, therefore, uncertainties associated with charging particles

# Naneum: Product Features

- Naneum products are state of the art and are designed for use in the field
- They are portable, robust and easy to use
- They operate untended and can be used as remote monitors.
- Naneum technology is based on Diffusion and Impaction Phenomena:
  - Reflects what happens in respiratory tract
  - Designed for ultra low pressure drop- no problems with volatiles and semi-volatiles
  - Does not involve charging- important in certain environments e.g. Carbon Nanotubes where issues such as multiple charging can distort distribution
  - No need for radiation sources or high voltage.
  - Affordable.